pure sense

CARE OF YOU

Dental C

WATER FLOSSER

PSWF500B

TABLE OF CONTENT

Safety precautions	4
Overview	7
Operation.....	8
Maintenance.....	11
Specifications.....	12

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SAFETY PRECAUTIONS

READ ALL OF THE INSTRUCTIONS BEFORE USE AND SAVE THEM FOR FUTURE REFERENCE.

WARNING! Risk of personal injury or product and/or property damage if you do not follow the precautions.

USAGE PRECAUTIONS

- Do not operate the device if you are unfamiliar with its functions or safety instructions. Otherwise, improper operation may cause injury.
- This product is designed for cleaning teeth and gums. Never use on eyes, nose, ears, or other body parts.
- This appliance is a personal care device and is not intended for use on multiple persons.
- If you are undergoing treatment for an oral care condition, consult your dental professional before use.
- If you have had oral surgery in the past 2 months, consult your dental professional before use.
- If you have a pacemaker or implanted medical device, contact your doctor or the manufacturer of the implanted device prior to use.
- This product is not suitable for children, infants, or individuals who are unable to operate it unless they are supervised or guided by an adult when using it.
- Keep this product out of the reach of children and infants, and children shall not play with the appliance.
- Cleaning and maintenance shall not be performed by children.
- Only use lukewarm or cool water in the water tank. Hot water may cause burns.
- Do not exert excessive pressure on the nozzle tip of the nozzle.
- Do not immerse this product in water to prevent malfunction, fire, electric shock, or explosion.
- Charge indoors only, away from heat sources (e.g., stoves, radiators).
- This appliance contains no user-serviceable parts. If damaged, stop

- use immediately and contact authorize support.
- Replace the nozzle every 6 months or sooner if worn /damage.
- Only use original accessories or consumables.
- Do not dismantle the device or modify its parts to avoid damaging the product, and the warranty will be invalidated.
- Any repairs to the product(s) must be done by qualified personnel. Do not try to repair the product(s) by yourself.

BATTERY PRECAUTIONS

- Disposal of a battery into fire or a hot oven, or mechanically crushing or cutting of a battery can result in an explosion.
- Leaving a battery in an extremely high temperature can result in an explosion or leakage of flammable liquid or gas; and a battery subjected to extremely low air pressure can result in an explosion or leakage of flammable liquid or gas.
- The unit may get hot during use, this is normal.
- This product shall not be thrown away or incinerated but disposed of in accordance with your local regulations concerning chemical wastes.
- A new battery with 100% capacity reduces its capacity by 20% after approximately 300 charging cycles.
- Disconnect the appliance from charging as soon as it is fully charged.
- Charge the battery every 3rd month if you do not use the product for a longer period of time. This is to ensure the battery performance.
- This product contains a lithium battery. When a device with a lithium battery is damaged (either once or several times repeatedly), the lithium-ion cells can be affected by a so-called 'thermal rush'. Thermal rush means that energy is released in the form of heat inside the battery and fire occurs quickly. It becomes a selfgenerating process that takes place explosively and is very difficult to extinguish.
- Do not apply excessive force on the product (do not hit or drop the product).
- Do not charge the product if it is very warm, to avoid charging a warm battery.
- Do not cover the product.



- Always inspect the product and stop using it if leakage, strange odors, color changes, deformation or any unusual occurrences are detected. If any leakage goes into your eyes, never rub them - clean your eyes with water and get medical help.
- Never disassembly or refit the unit. Repairs to the equipment shall only be performed by qualified personnel. Improper repairs may place the user at serious risk.

SYMBOL	EXPLANATION
	Read all of the instructions before use and save them for future reference.
	This manual does not provide any warranty or related declarations. All information in this manual including the product specification, pictures and literatures is just for reference. NetOnNet AB reserves the rights to update contents of this manual without further notice.
	This manual is the copyright of NetOnNet AB. No part of this manual may be revised, copied or transmitted in any way without written permission from NetOnNet AB.
	The CE mark indicates that this product has been manufactured according to current regulations for electromagnetic fields, which means that it is safe to use if used according to safety instructions. The product is manufactured according to current regulations within the EU legislation, limiting the use of hazardous substances in electrical and electronic equipment.
	This product and the accessories packed together are the applicable product to the WEEE directive except batteries. Please recycle any materials in accordance with your local recycling regulations. Batteries should never be thrown away or incinerated but disposed of in accordance with your local regulations concerning chemical wastes.
	This product's packaging materials are recyclable and can be reused. Please recycle any materials in accordance with your local recycling regulations.
IPX7	Waterproof
	This symbol identify equipment meeting the safety requirements specified for Class II equipment according to IEC 61140.
	The exclamation point within an equilateral triangle is intended to alert the user to the presence of important instructions. If these instructions are not followed, there is risk of personal injury or product and/or property damage.

OVERVIEW

PARTS

1. Nozzle
2. Nozzle release button
3. Power switch
4. Mode switch
5. Pulse mode indicator light
6. Pulse button
7. Power indicator
8. Strong mode
9. Normal mode
10. Soft mode
11. Pulse mode
12. Water tank
13. Charging port

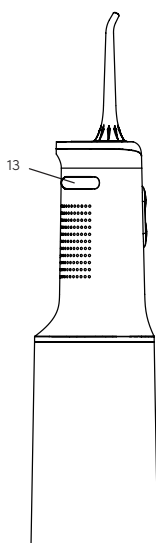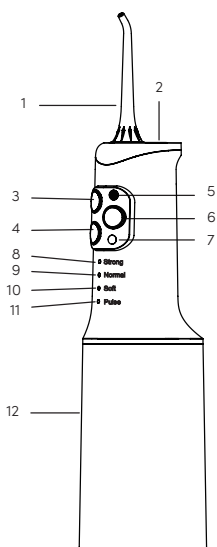
OPERATION

CHARGE THE WATER FLOSSER


1. During charging, place the main unit on a horizontal surface in a reclined position to avoid the risk of the product tipping over.
2. Make sure there is no liquid left in the water tank before charging.
3. Use the included USB-C cable for charging. The charging time is approximately 6 hours for 0-100%.
4. When the charging cable is plugged in, the green light will breathe and blink (1 second/time) until the battery is fully charged, at which time the green light will remain on.
5. During operation, the power indicator is green when the battery level is above 10% (normal operation) and red when the battery level is below 10% (low battery).

Note: If the battery is normal, the light will turn off when the unit is turned off; if the battery is low, the red light will flash after the unit is turned off and then turn off after 5 seconds.

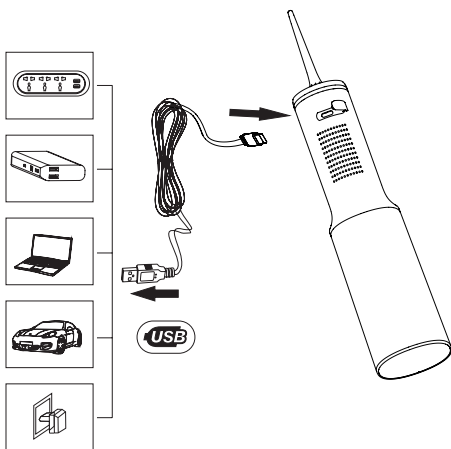
GET STARTED:


1. Open the water tank lid, fill the tank with water not exceeding 40°C, and close the lid securely.

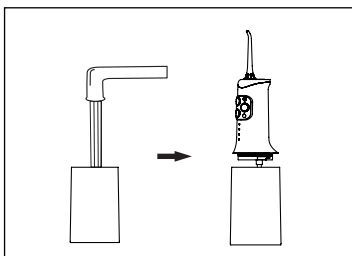
2. For the first use, add a small amount of clean water to the outlet.
3. Insert the nozzle into the nozzle port on the main unit until a clicking sound is heard, indicating that it is fully inserted.

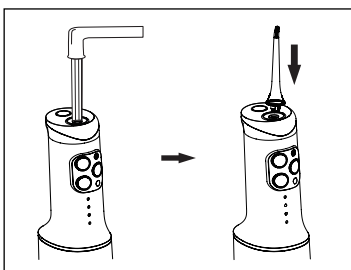4. Installation of the water tank
Rotate the main unit to align it with the corresponding slot on the water tank, then push the main unit forward to lock it in place.

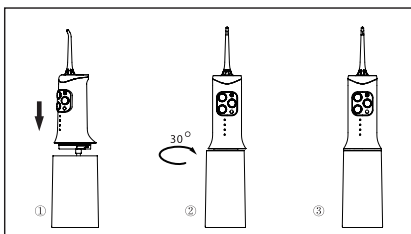5. Removal of the water tank
Rotate the main unit backward while keeping the water tank in position, misalign the main unit and tank, and take out the main unit.

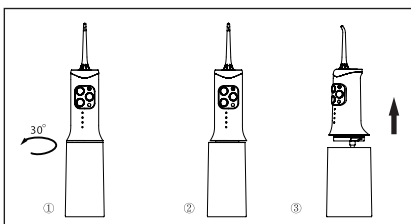
USAGE:


1. In the power-off state, press the mode button once will switch the gear according to the sequence of Strong → Normal → Soft → Pulse.
2. Mode description:

Mode	Suitable use
Strong	Deep cleaning of food debris stuck between teeth.
Normal	Individuals with sensitive gums and daily cleaning.
Soft	Individuals with sensitive gums and first-time users.
Pulse	Alternating water flow strength with pulsed massage for gums.

3. Press the power switch button once will start the motor to work, and pressing it again will stop the motor. Each gear has a 2-minutes timer.
4. You can also press and hold pulse button to spray water, water spray stops when releasing the button.

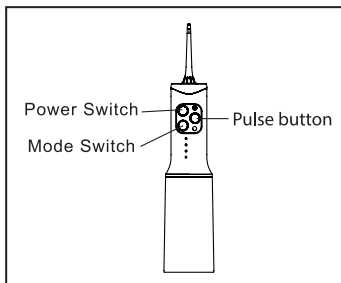5. Keep the nozzle tip vertically aligned with your teeth. Slightly open your mouth to allow the water to flow out. Clean your teeth.

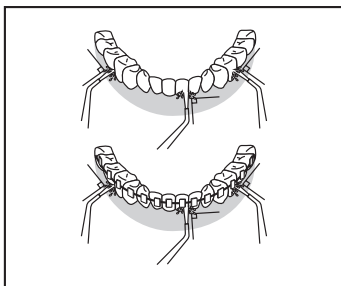
6. Spray water onto your teeth and gums at a 90-degree angle.
7. Guide the water to slowly flow over your teeth.
8. Spray water onto the gum tissue along the gum line with the nozzle tip at a 90-degree angle.
9. Thoroughly clean your teeth and the areas around braces, crowns, and bridges.


MAINTENANCE

CLEANING:

1. Press the power switch button to turn off the device.
2. Press the nozzle lock button, and then remove the nozzle by lifting it up.

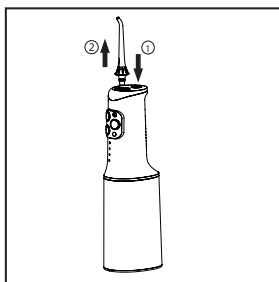3. Remove the water tank and turn it upside down to discard the water inside.

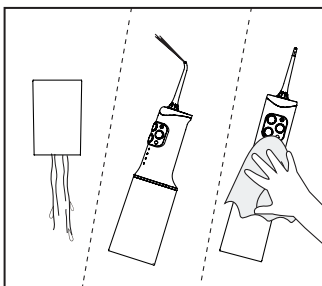4. Press the power switch button to turn on the device and allow the water tube to fully drain.
5. Wipe the device with a dry cloth.

REMOVING HARD WATER DEPOSITS / GENERAL CLEANING

- Hard water deposits may accumulate inside the device depending on the mineral content of the water. If ignored, this can affect performance.
1. To clean internal components: fill the water tank with warm water and add 1 tablespoon of white vinegar.
 2. Direct the handle and nozzle towards the sink. Turn on the device and let it run until the tank is empty.
 3. Then, refill the tank with clean warm water and repeat the process to flush it out.
 4. Repeat this operation every 1 to 3 months to ensure optimal performance.

STORAGE

- The product should be stored away from dust, dirt and water sources such as tubs and sinks, or in a travel case.
- Store the product out of reach of children.

SPECIFICATIONS

Model name	PSWF500B
Battery	3.7V, 2000mAh, 7.4Wh
Charging time	6 hours
Working time	40 minutes (Normal mode)
Water tank capacity	300 ml
Spray frequency	1600 times/minute
IP class	IPX7
Cable length	95 cm
Size	71 x 61 x 231 mm
Weight	295 g
Included in the package	- Water flosser - USB-A to USB-C charging cable - Nozzles x4 - User manual

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